

Reference 620

Potentiostat/Galvanostat/ZRA

- + Fast CV
- + EIS To 5 MHz
- + Ultra Low Current
- + Floating Operation
- + Low Noise



GAMRY
INSTRUMENTS

Reference 620

Potentiostat/Galvanostat/ZRA




The **Reference 620** has been designed to bring you maximum performance for demanding applications. It is ideally suited for:

- + Corrosion Measurement
- + Paints & Coatings
- + Bioelectrochemistry & Sensor Development
- + Physical & Analytical Electrochemistry
- + Energy Devices

+ High Performance EIS

EIS from 10 μ Hz to 5 MHz. Our specially designed cables reduce mutual inductance increasing the low impedance bandwidth. Techniques include single sine Potentiostatic, Galvanostatic, and Hybrid EIS. Also included are multi-sine techniques for Potentiostatic and Galvanostatic. Our unique power-leveling algorithm improves signal-to-noise ratios and reduces acquisition time.

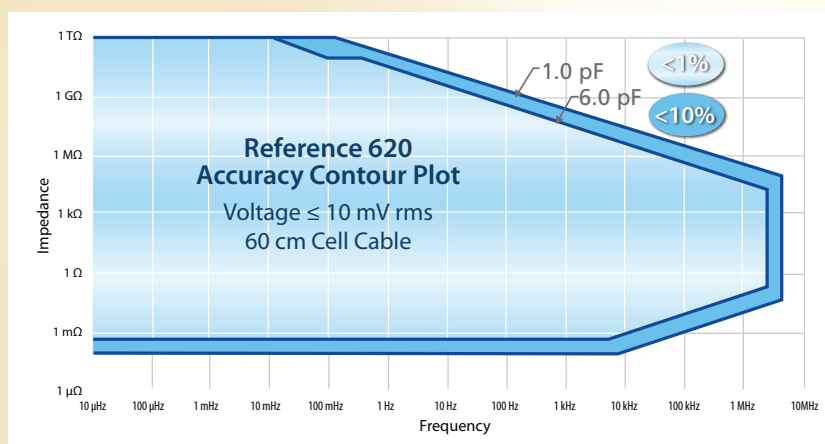


+ Ultra-Low Noise

The Reference 620 enjoys an intrinsic noise level of <2 μ V rms, thanks to a well-designed electronic layout, components selected to balance speed and noise, and intelligent analog filtering and shielding. The instrument can oversample the signals, then average to smooth a slowly changing, but noisy, signal.

+ Digital Signal Processing

All Gamry potentiostats employ digital signal processing (DSP) technology, allowing for oversampling and averaging in order to improve signal-to-noise ratios and provide accurate capacitance measurements. Our instruments have three sampling modes – Fast, Noise Reject, and Surface. Fast corresponds to sampling at the end of each step. Noise reject oversamples and averages during the last 20% of a step. Surface mode oversamples and averages during the entire step thereby ensuring no lost charge for an accurate capacitance measurement.

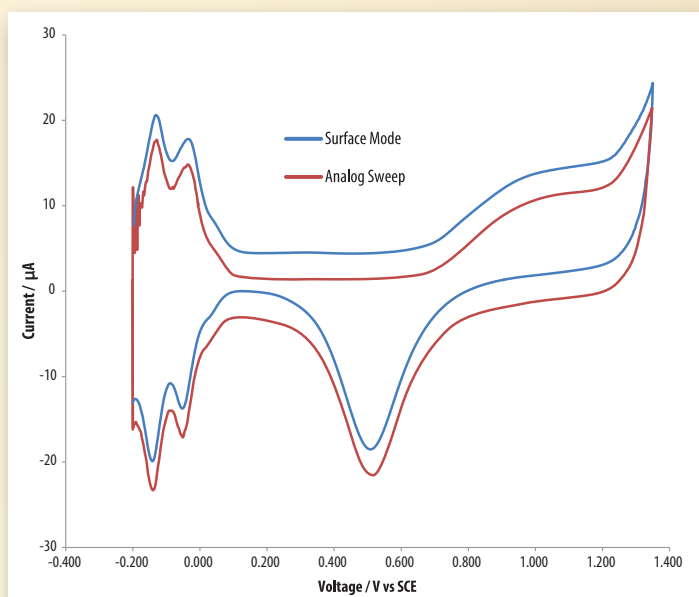

Comparison of Surface Mode vs Analog Sweep. Pt WE in 1 M H₂SO₄. Surface mode current offset 3 μA for easy comparison.

+ Floating Operation

Gamry provides the highest electrical isolation allowing for measurements on grounded cells and electrodes. Electrical isolation also allows for coupling to other instruments such as TEMs and SECMs.

+ Small Footprint

Only 9x19x27 (WxHxD) cm and 3 kg.

+ Current Interrupt & Positive Feedback iR Compensation

The Reference 620 has both current interrupt and positive feedback modes of iR Compensation.

+ Optional Inputs & Outputs

The Reference 620 gives you several input and output options. An additional voltage can be read via a BNC input. An external signal can be fed into the control amplifier, allowing you to input external signals right into your cell. A 15-pin connector is provided to allow for digital inputs and outputs for TTL level triggering of or with external devices. This connector also provides a voltage out, allowing control of external devices such as a rotating electrode setup.

Reference 620

Potentiostat/Galvanostat/ZRA

SPECIFICATIONS

Potentiostat	Yes
Galvanostat	Yes
Zero Resistance Ammeter	Yes
Cell Connections	2, 3, 4, 5, or 6
Isolated from earth	Yes

SYSTEM

Max. Current	± 600 mA
Current Ranges	11 (60 pA-600 mA), 13 with 10X and 100X gain
Max. Applied Potential	± 12 V
Rise Time	< 250 ns
Noise and Ripple	< 2 μ V rms (typical)
Min. Time Base	3.333 μ s
Min. Potential Step	12.5 μ V

EIS MEASUREMENT

Frequency Range	10 μ Hz – 5 MHz
Impedance Accuracy	See accuracy contour map
Max AC Amplitude	3 V max 600 mA max

CONTROL AMP

Compliance Voltage	$> \pm 22$ V
Output Current	$> \pm 600$ mA
Speed Settings	5

ELECTROMETER

Input Impedance	$> 10^{14}$ Ω
Input Current	< 10 pA
Bandwidth (-3dB) (typical)	> 15 MHz
Common Mode Rejection Ratio (CMRR)	> 65 dB (1 MHz)

